

Model No. IMRX44680

QUESTIONS?

As a manufacturer, we are committed to providing complete customer satisfaction. If you have questions, or if there are missing parts, we will guarantee complete satisfaction through direct assistance from our factory.

TO AVOID UNNECESSARY DELAYS, PLEASE CALL DIRECT TO OUR TOLL-FREE CUSTOMER HOT LINE. The trained technicians on our customer hot line will provide immediate assistance, free of charge to you.

CUSTOMER HOT LINE:

1-800-999-3756

Mon.–Fri., 6 a.m.–6 p.m. MST

CAUTION

Read all precautions and instructions in this manual before using this equipment. Save this manual for future reference.

IMAGE[®]

RENEW 90/90

Back Therapy Kit



Patent Pending

USER'S MANUAL

Before You Begin

Thank you for selecting the innovative IMAGE RENEW 90/90 Back Therapy Kit. The 90/90 Back Therapy Kit is designed to reduce stress on your back by creating a 90° angle between your torso and your thighs. It also provides heat and vibration massage. At the same time, the 90/90 is portable and so versatile that it can be used at home, in the office or in your car.

For your benefit, read this manual carefully before using the 90/90. If you have additional questions, please call our Customer Service Department at the number shown to the left. To help us assist you, please note the product model number before calling. The model number is IMRX44680.

Important Safety Precautions

WARNING: To reduce the risk of burns, fire, electric shock or serious injury, read these precautions and all instructions in this manual before using the 90/90.

1. Read all precautions and instructions in this manual before using the 90/90.
2. It is the responsibility of the owner to ensure that all users of the 90/90 are adequately informed of all precautions.
3. Always turn off the hand control and unplug the power cord immediately after using and before cleaning the massage pad.
4. Never leave the massage pad unattended while it is plugged in.
5. Do not use the massage pad under a blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to persons.
6. Close supervision is necessary when the massage pad is used by, on or near invalids or disabled persons.
7. Keep children under the age of 12 and pets away from the 90/90 at all times.
8. Use the 90/90 only for its intended use as described in this manual.
9. Carefully examine the massage pad before each use. Never operate the massage pad if it is not working properly, if it is damaged or if the power cord or plug is damaged. If any of these conditions exists, call the toll-free telephone number on the front of this manual.
10. Keep the power cord away from heated surfaces.
11. The Back Therapy Kit is intended for home use only. Do not use the 90/90 in a commercial, rental or institutional setting.

DANGER


12. Keep pins and other sharp objects away from the massage pad. Do not crush the massage pad. Avoid sharp folds in the massage pad.
13. Use the massage pad indoors. Never use the massage pad outdoors. Keep the massage pad dry, away from water and moisture.
14. The vibration motor and the heating element do not shut off automatically. Temperatures sufficiently high to cause burns may occur. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the massage pad frequently to reduce the risk of blistering. Limit use to 20 minutes at a time to avoid bruises to the skin.

WARNING: Read all instructions before using. ICON assumes no responsibility for personal injury or property damage sustained by or through the use of this product.

Using the 90/90 Back Therapy Kit

Orientation Diagram

This drawing shows the major parts of the Back Therapy Kit.

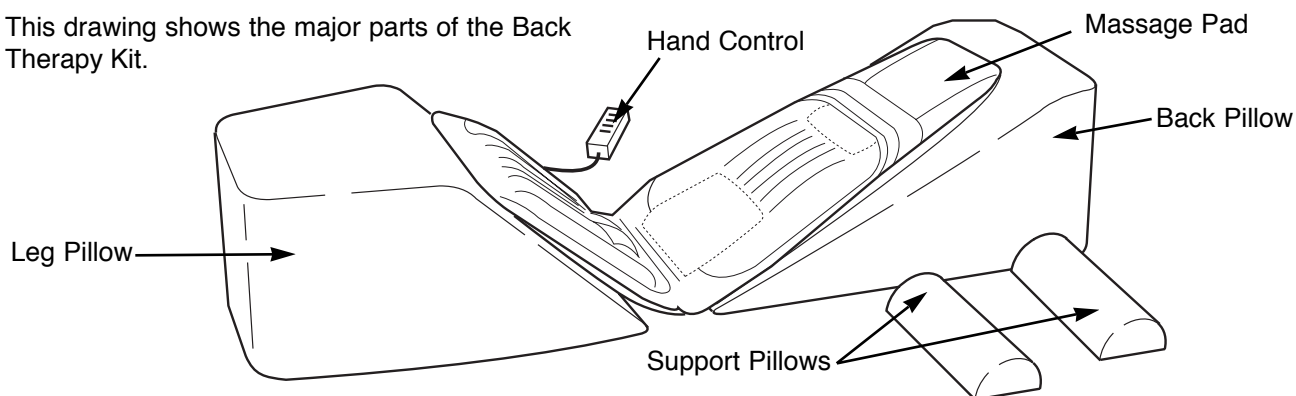

Using all Parts of the Back Therapy Kit



Place the back pillow and the leg pillow on the floor with the massage pad on top. Plug in the massage pad as described on the next page. This combination gives you maximum support and pain relief by creating a 90° angle between torso and thighs.

Using the Back Pillow and Massage Pad



Place the back pillow on the floor with the massage pad on top. This combination supports your back and neck while allowing your legs to be stretched out comfortably.

Using the Leg Pillow and Massage Pad



Place the leg pillow on the floor with the massage pad on top. The massage pad can be under your back as in the previous combinations or under your legs only. This combination raises your legs above the level of your heart, thereby increasing circulation.

Using the Massage Pad Only



The massage pad can be used alone either on the floor or on a bed. **CAUTION:**

Note: The two support pillows can be placed under your neck or knees for any of the combinations described above.

Using the Massage Pad on a Chair



The massage pad can be used on most chairs, such as a recliner or office chair. Strap the massage pad onto the chair with the elastic straps on the back of the pad.

By using the included auto adapter, the pad can be used in your car to relieve fatigue and discomfort during a long drive.

Using the Hand Control

The hand control attached to the massage pad must be plugged into a power source for the pad to function. The hand control can be plugged into a wall outlet with the included AC adapter or into the cigarette lighter of a car with the included auto adapter. When it is not being used, the hand control can be stored in the pocket on the massage pad.

The hand control has the following buttons:

Power On/Off. Turns the power on and off. A red light is visible when the power is on.


Heat On/Off. Turns the heating element on and off. A red light is visible when the heat is on. The heating element has a built-in thermostat to prevent overheating.

Zone 1. Controls the vibrations in the neck area of the massage pad. Note: All of the zone buttons have a low, medium and high setting, allowing you to regulate the intensity of the massage.

Zone 2. Controls the vibrations in the upper back area.

Zone 3. Controls the vibrations in lower back area.

Zone 4. Controls the vibrations in the seat area.

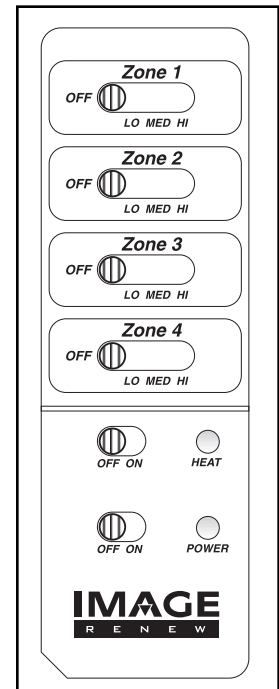

Ordering Replacement Parts

To order replacement parts, call our Customer Service Department toll-free at 1-800-999-3756, Monday through Friday, 6 a.m. until 6 p.m. Mountain Time (excluding holidays). To help us assist you, please be prepared to give the following information:

- The MODEL NUMBER of the product (IMRX44680).
- The NAME of the product (Image Renew 90/90 Back Therapy Kit).

Limited Warranty

ICON Health & Fitness, Inc. (ICON), warrants this product to be free from defects in workmanship and material, under normal use and service conditions, for a period of ninety (90) days from the date of purchase. This warranty extends only to the original purchaser. ICON's obligation under this warranty is limited to replacing or repairing, at ICON's option, the product at one of its authorized service centers. All products for which warranty claim is made must be received by ICON at one of its authorized service centers with all freight and other transportation charges prepaid, accompanied by sufficient proof of purchase. All returns must be pre-authorized by ICON. This warranty does not extend to any product or damage to a product caused by or attributable to freight damage, abuse, misuse, improper or abnormal usage or repairs not provided by an ICON authorized service center, products used for commercial or rental purposes, or products used as store display models. No other warranty beyond that specifically set forth above is authorized by ICON.

ICON is not responsible or liable for indirect, special or consequential damages arising out of or in connection with the use or performance of the product or damages with respect to any economic loss, loss of property, loss of revenues or profits, loss of enjoyment or use, costs of removal, installation or other consequential damages of whatsoever nature. Some states do not allow the exclusion or limitation of incidental or consequential damages. Accordingly, the above limitation may not apply to you.

The warranty extended hereunder is in lieu of any and all other warranties and any implied warranties of merchantability or fitness for a particular purpose is limited in its scope and duration to the terms set forth herein. Some states do not allow limitations on how long an implied warranty lasts. Accordingly, the above limitation may not apply to you.

This warranty gives you specific legal rights. You may also have other rights which vary from state to state.

ICON HEALTH & FITNESS, INC., 1500 S. 1000 W., LOGAN, UT 84321-9813