

Limited Warranty

ICON Health & Fitness, Inc. (ICON), warrants this product to be free from defects in workmanship and material, under normal use and service conditions, for a period of ninety (90) days from the date of purchase. This warranty extends only to the original purchaser. ICON's obligation under this warranty is limited to replacing or repairing, at ICON's option, the product at one of its authorized service centers. If the product fails to perform when used properly because of a manufacturing defect, simply return the product to the retail location where it was purchased, or send it postage prepaid to the address below within 90 days following the original date of purchase. This warranty does not extend to any product or damage to a product caused by or attributable to freight damage, abuse, misuse, improper or abnormal usage or repairs not provided by an ICON authorized service center, products used for commercial or rental purposes, or products used as store display models. No other warranty beyond that specifically set forth above is authorized by ICON.

ICON is not responsible or liable for indirect, special or consequential damages arising out of or in connection with the use or performance of the product or damages with respect to any economic loss, loss of property, loss of revenues or profits, loss of enjoyment or use, costs of removal, installation or other consequential damages of whatsoever nature. Some states do not allow the exclusion or limitation of incidental or consequential damages. Accordingly, the above limitation may not apply to you. The warranty extended hereunder is in lieu of any and all other warranties and any implied warranties of merchantability or fitness for a particular purpose is limited in its scope and duration to the terms set forth herein. Some states do not allow limitations on how long an implied warranty lasts. Accordingly, the above limitation may not apply to you.

This warranty gives you specific legal rights. You may also have other rights which vary from state to state.

ICON Health & Fitness, Inc., 1500 S. 1000 W., Logan, UT 84321

Part No.161143 R1199A Printed in USA © 1999 ICON Health & Fitness, Inc.

Limited Warranty

ICON Health & Fitness, Inc. (ICON), warrants this product to be free from defects in workmanship and material, under normal use and service conditions, for a period of ninety (90) days from the date of purchase. This warranty extends only to the original purchaser. ICON's obligation under this warranty is limited to replacing or repairing, at ICON's option, the product at one of its authorized service centers. If the product fails to perform when used properly because of a manufacturing defect, simply return the product to the retail location where it was purchased, or send it postage prepaid to the address below within 90 days following the original date of purchase. This warranty does not extend to any product or damage to a product caused by or attributable to freight damage, abuse, misuse, improper or abnormal usage or repairs not provided by an ICON authorized service center, products used for commercial or rental purposes, or products used as store display models. No other warranty beyond that specifically set forth above is authorized by ICON.

ICON is not responsible or liable for indirect, special or consequential damages arising out of or in connection with the use or performance of the product or damages with respect to any economic loss, loss of property, loss of revenues or profits, loss of enjoyment or use, costs of removal, installation or other consequential damages of whatsoever nature. Some states do not allow the exclusion or limitation of incidental or consequential damages. Accordingly, the above limitation may not apply to you. The warranty extended hereunder is in lieu of any and all other warranties and any implied warranties of merchantability or fitness for a particular purpose is limited in its scope and duration to the terms set forth herein. Some states do not allow limitations on how long an implied warranty lasts. Accordingly, the above limitation may not apply to you.

This warranty gives you specific legal rights. You may also have other rights which vary from state to state.


ICON Health & Fitness, Inc., 1500 S. 1000 W., Logan, UT 84321

Part No.161143 R1199A Printed in USA © 1999 ICON Health & Fitness, Inc.

Hot&Cold TherapyPack

with Patented
COMFORTEMP

WLRX12590
WLRX25590

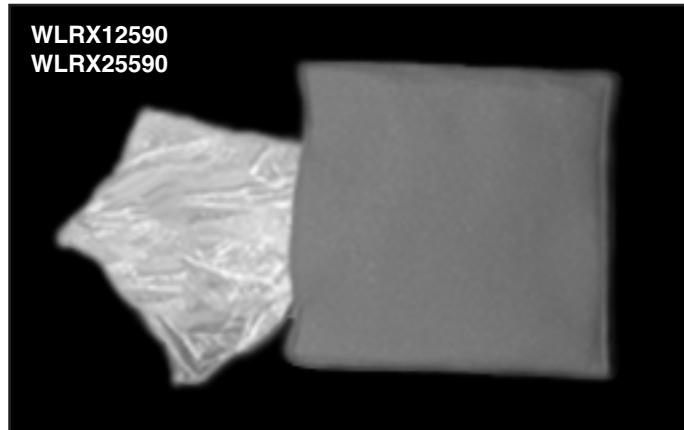


Thank you for purchasing the WESLO® Hot & Cold Therapy Pack. For your benefit, please read this User's Guide carefully before using this product.

WESLO is a registered trademark of
ICON Health & Fitness, Inc.

Hot&Cold TherapyPack

with Patented
COMFORTEMP

WLRX12590
WLRX25590

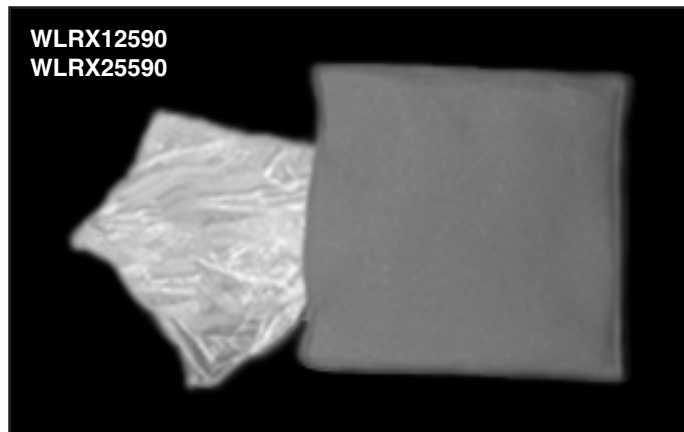

Thank you for purchasing the WESLO® Hot & Cold Therapy Pack. For your benefit, please read this User's Guide carefully before using this product.

WESLO is a registered trademark of
ICON Health & Fitness, Inc.

Preparing the Gel Pack for Cold or Hot Therapy

Directions for Cold Therapy



1. Place the gel pack in the freezer for a minimum of four hours.
2. Remove the gel pack from the freezer. Insert the gel pack into the Hot/Cold Pillow Case. The blue side of the Pillow Case denotes cold, and should be placed against your body when using cold therapy.
3. Apply the cold therapy pack to the desired area of your body.

Directions for Hot Therapy




1. Using either the stove or microwave oven, boil enough water to entirely cover the gel pack.
2. Place the gel pack into an appropriate sized container. Pour the heated water over the gel pack so that it is completely submerged for approximately four minutes.
3. Remove the gel pack from the heated water. **Caution: The water and gel pack are hot. Be careful not to burn yourself while removing the gel pack from the water.** Insert the gel pack into the Pillow Case. The red side of the Pillow Case denotes heat, and should be placed against your body.

4. Apply the hot therapy pack to the desired area of your body.


Washing the Pillow Case

To wash the Hot & Cold Therapy Pack pillow case, first remove the gel pack, as shown. Wash in warm water and tumble dry. **Do not wash the gel pack.**

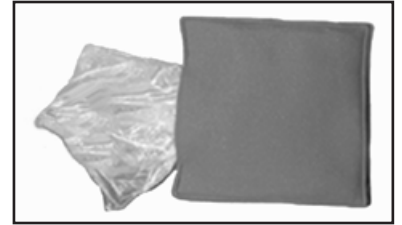

ComforTemp® DCC™ prolongs Hot/Cold Therapy

Hot & Cold Therapy Pack features ComforTemp® foam. ComforTemp® DCC™ (Dynamic Climate Control) contains Thermasorb® additives that absorb and release heat to retain hot or cold temperatures longer than standard gel packs. Prolong the benefits of hot/cold therapy with ComforTemp® foam.

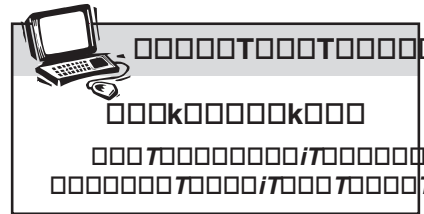

Preparing the Gel Pack for Cold or Hot Therapy

Directions for Cold Therapy



1. Place the gel pack in the freezer for a minimum of four hours.
2. Remove the gel pack from the freezer. Insert the gel pack into the Hot/Cold Pillow Case. The blue side of the Pillow Case denotes cold, and should be placed against your body when using cold therapy.
3. Apply the cold therapy pack to the desired area of your body.

Directions for Hot Therapy



1. Using either the stove or microwave oven, boil enough water to entirely cover the gel pack.
2. Place the gel pack into an appropriate sized container. Pour the heated water over the gel pack so that it is completely submerged for approximately four minutes.
3. Remove the gel pack from the heated water. **Caution: The water and gel pack are hot. Be careful not to burn yourself while removing the gel pack from the water.** Insert the gel pack into the Pillow Case. The red side of the Pillow Case denotes heat, and should be placed against your body.

4. Apply the hot therapy pack to the desired area of your body.

Washing the Pillow Case

To wash the Hot & Cold Therapy Pack pillow case, first remove the gel pack, as shown. Wash in warm water and tumble dry. **Do not wash the gel pack.**



ComforTemp® DCC™ prolongs Hot/Cold Therapy

Hot & Cold Therapy Pack features ComforTemp® foam. ComforTemp® DCC™ (Dynamic Climate Control) contains Thermasorb® additives that absorb and release heat to retain hot or cold temperatures longer than standard gel packs. Prolong the benefits of hot/cold therapy with ComforTemp® foam.

